

INJECT NEW VITALITY INTO THE SKIN





The world's first annular "micro- pulse ultrasound" collaborative technology, ultrasonic linear lattice focusing, can maximize the operation effect; In addition to lifting, it also brings extra firmness Compared with traditional ultrasound mode, multiple modes can effectively reduce pain



Multidimensional flat handle probe

Facial lift and skin firmness. Display the current handle with which probe handle. When the handle is lifted, it automatically switches the handle operation interface.



Multidimensional flat handle probe


Reshape the probe. Display the current handle with which probe handle. When the handle is lifted, it automatically switches the handle operation interface.



360 ° annular handle probe

Facial lift and skin firmness. Display the current handle with which probe handle. When the handle is lifted, it automatically switches the handle operation interface (Figure 3).

When you pick up the first handle, the instrument automatically recognizes and converts the corresponding handle probe operation interface. Whether you pick up another handle or two handles later, the instrument identifies and displays the operation interface for picking up the first handle. If you need to change the probe operation interface, first place all handles on the instrument, and then pick up the handle you need to operate, The interface automatically recognizes and converts the probe corresponding to the operation interface when you pick up the first handle.

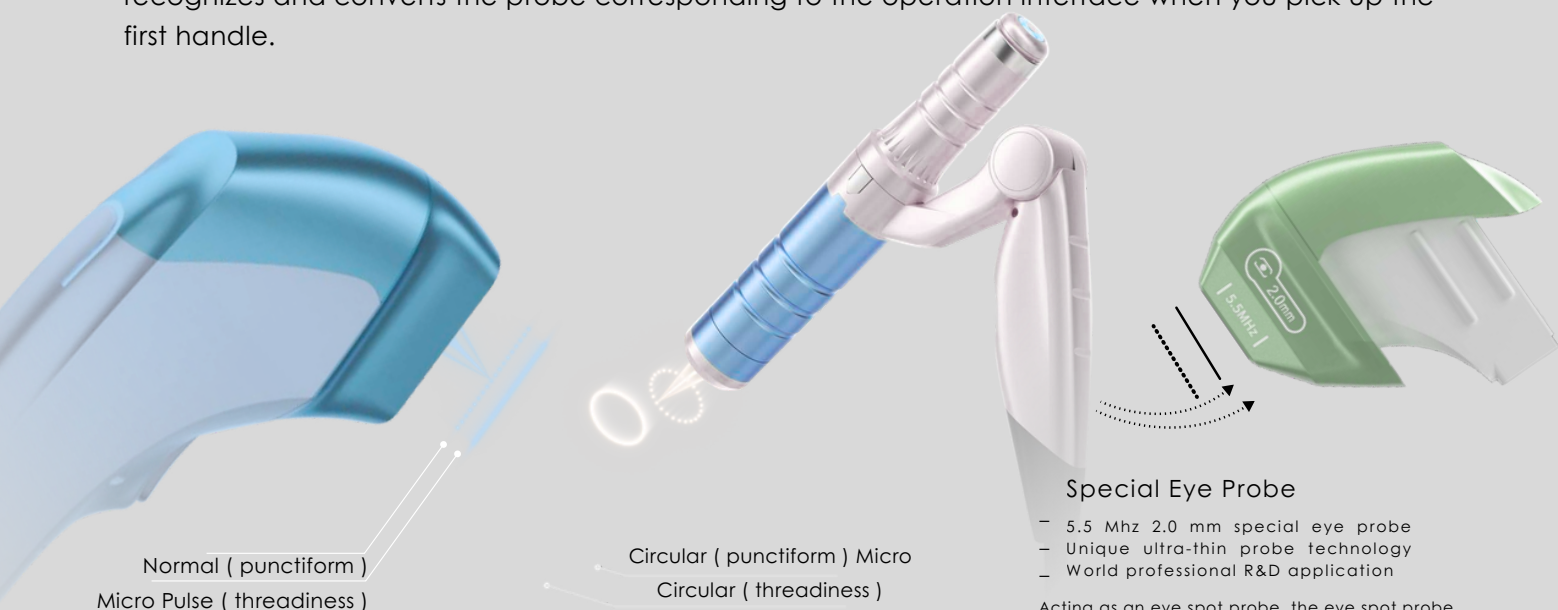

Special Eye Probe

- 5.5 Mhz 2.0 mm special eye probe
- Unique ultra-thin probe technology
- World professional R&D application

Acting as an eye spot probe, the eye spot probe refines and solves problems such as shallow fine lines, deep genuine lines, dark circles, and eye bags.

TECHNICAL SPECIFICATIONS

Parameter	Specification
Product Name	13D HIFU Cartridges
Handle Cartridges	3 pcs (included)
Length	5mm / 10mm / 15mm / 20mm / 25mm (adjustable)
Cartridge Shots	10,000 shots per cartridge
Power	650W (High-Intensity Focused Ultrasound)
Frequency	2MHz / 4MHz / 5.5MHz (adjustable)
Modes	Normal Mode / Multi-Pulse (MP) Mode

The Body Probe Is Applied To The Skin

Normal (Punctate)

.....

Micro Pulse (Linear)

—————

Apply gel to the skin before operation
Operating handle vertical skin DO not
leave the skin during operation

T-6.0 (body) 2MHz


Apply to larger body parts like arms or thighs
to tighten skin and improve the appearance
of cellulite.

T-9.0 (body) 2MHz

Specifically targets loose abdominal skin and
improves body contours.

T-13 (body) 2MHz

For treatment of larger areas: such as
the flanks, back, and thighs where fat
layers are thick

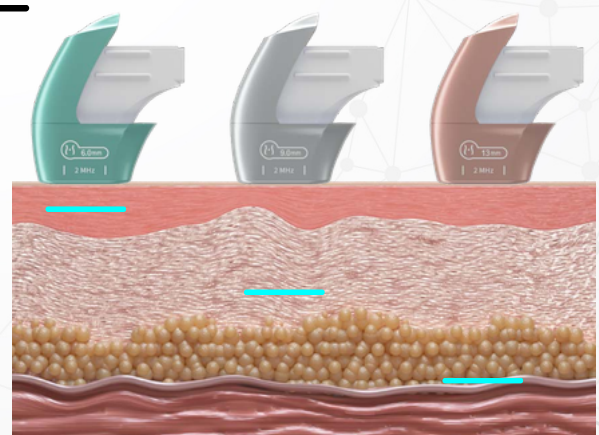

Achieve comprehensive slimming and shaping.

